

RP105 antibody [MHR73] (FITC)

Cat. No. GTX75888

宿主	Mouse
克隆	Monoclonal
同种型	IgG1
实验应用	FCM
种属反应	Human

实验应用

应用说明

*最佳稀释倍数与浓度应由研究人员确认

Suggested dilution	Recommended dilution
FCM	Neat-1/10

Note : CD180 is expressed at the cell surface in association with MD-1 and it may also be found free within the cell. Use 10 μ l of the suggested working dilution to label 10⁶ cells in 100 μ l.

以下为常规应用缩写的中文注解

WB: 免疫印迹

ICC/IF: 细胞染色

IHC-P: 石蜡切片

IHC-Fr: 冰冻切片


属性

形式	Liquid
存储溶液	PBS, 1% BSA
保存剂	0.09% Sodium azide
存放说明	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles. Protect from light.
浓度	0.1 mg/ml (Please refer to the vial label for the specific concentration.)
偶联	Fluorescein isothiocyanate (FITC) Wavelength
RRID	AB_377097
注意事项	仅供实验室使用。不适用于人类或动物的任何临床、治疗或诊断用途。不适合动物或人类食用。



For full product information, images and publications, please visit our [website](#).

產品圖片

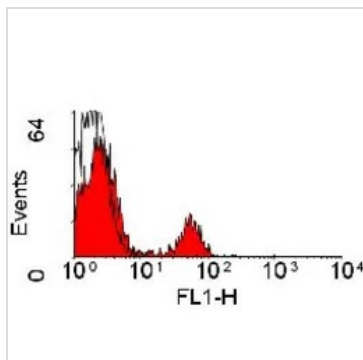

GTX75888 FCM Image

FACS analysis of human peripheral blood lymphocytes using GTX75888 RP105 antibody [MHR73] (FITC).



For full product information, images and publications, please visit our [website](#).